

# Laryngeal Mask Airways

## DESCRIPTION

m|devices Laryngeal Mask Airways are designed to effectively protect the air ways of patients with spontaneous breathing or positive pressure breathing. They are made from 100% medical grade silicone and its smooth one-piece design with guide markings on the handle, assist with easy insertion and positioning. The translucent tubing allows clear visibility and the soft seal cuff, which is permeable to nitrous oxide, reduces pressure rise and minimises potential trauma. Our masks feature a 15mm connector and a pilot balloon to control the cuff pressure. m|devices Laryngeal masks are also available in 7 colour coded sizes for easy and fast identification and are stored in single, sterile packaging.

## FEATURES & BENEFITS

- Pilot Balloon with non return valve
- 100% Medical grade silicone
- Smooth one-piece design
- Guide markings
- Colour coded for easy size selection
- Translucent tubing
- Soft seal cuff
- Peel open pouch
- Single patient use only


## PRODUCT DATA

- Maximum Cuff Inflation guide -  
Size 1. - 4ml, Size 1.5- 7ml,  
Size 2-10ml, Size 2.5-14ml,  
Size 3- 20ml, Size 4- 30ml,  
Size 5- 40ml
- Not made with Natural Rubber Latex
- Not made with DEHP
- Sterile EO
- Material: Silicone
- Shelf Life: 3 Years
- ARTG Number: 210988

## STANDARDS

- EN ISO 14155-1:2003 Clinical investigation for medical devices for human
- EN 866:1999 Biology system for testing sterilizer and sterilisation process
- EN 1174:1996 Sterilisation of medical Devices

## PRODUCT DIAGRAM

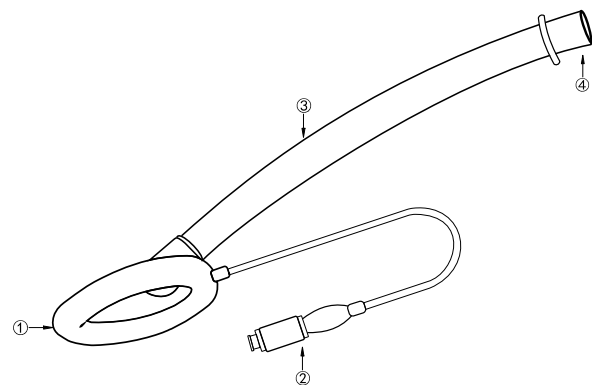

ITEM	DESCRIPTION
1	Cuff - Colour Coded
2	Pilot Balloon
3	Standard Tube
4	Connector

# Laryngeal Mask Airways

Laryngeal Mask Airways			
RE-ORDER	DESCRIPTION	SIZE	CARTON/QTY
AN030000	Laryngeal Mask Airway	1	10
AN030001	Laryngeal Mask Airway	1.5	10
AN030002	Laryngeal Mask Airway	2	10
AN030003	Laryngeal Mask Airway	2.5	10
AN030004	Laryngeal Mask Airway	3	10
AN030005	Laryngeal Mask Airway	4	10
AN030006	Laryngeal Mask Airway	5	10