

# Cervical Dilating Balloon

## DESCRIPTION

m|devices Cervical Dilating Balloon and Universal Securement Device has been specifically designed to be used as a mechanical dilation device for cervical ripening. The round closed distal end meets infection control principles by maintaining a closed system and clear markings for balloon capacity on the balloon lumen and printed batch have all been considered to assist the clinician with ease of use. The Universal Securement Device (optional) is included to apply light traction post balloon insertion.

## FEATURES & BENEFITS

- Silicone material assists with an easy insertion and does not require the aid of an introducer
- Balloon lumen is distinguished by a yellow coloured valve and is clearly labeled with the maximum volume capacity.
- High contrast position markings
- Round closed end maintains infection control principles which can often be a concern with a Foley catheter, when used in this application
- Batch number printed for ease of identification


## PRODUCT DATA

- 80mL balloon capacity
- Radiopaque
- Size: 20Fr, Length 40cm
- Not Made with Natural Rubber Latex
- Not made with DEHP
- Materials: Balloon Catheter: Silicone
- ARTG 307304
- Shelf life: 3 years

## STANDARDS

- ISO standards- ISO 11135:2014- Sterilisation of health care products. Requirements for the development , validation and routine control of a sterilization process for medical devices.
- ISO 14971:2012 Medical devices-Application of risk management to medical devices

## PRODUCT DIAGRAM

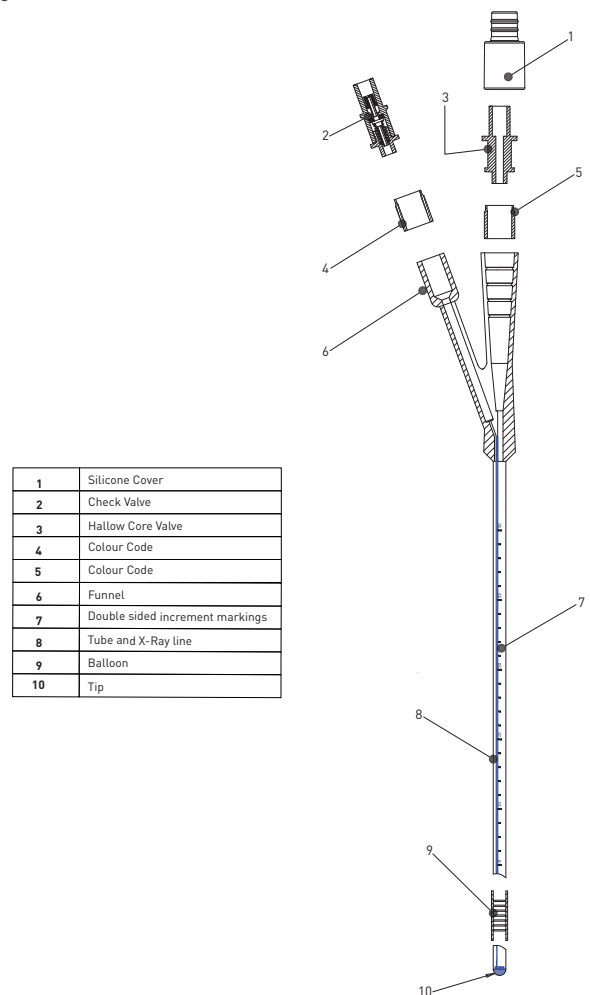

# Cervical Dilating Balloon

Cervical Dilating Balloon			
RE-ORDER	DESCRIPTION	SIZE	CARTON/QTY
<b>OB010020</b>	Cervical Dilating Balloon	80mL balloon	10