

Cervical Dilating Balloon Introducer

DESCRIPTION

m|devices Cervical Dilating Balloon Introducer has been designed with clinicians to assist with the insertion of the m|devices Cervical Dilating Balloon. The malleable nylon material allows the clinician to manipulate a structured shape for a simple and efficient insertion of the Cervical Dilating Balloon. The J-loop handle allows the clinician to secure the Introducer within the Cervical Dilating Balloon once in place to minimise movement. Removal of the introducer is to simply reverse bend the J-loop handle and the Introducer can be quickly removed from the m|devices Cervical Dilating Balloon.

FEATURES & BENEFITS

- Malleable material to meet user requirements for insertion
- J-loop handle enables the introducer to be securely in place during use
- Can be manipulated into any angle/shape as needed
- Easy to insert and remove from the m|devices Cervical Dilating Balloon


PRODUCT DATA

- Material: Blue Nylon
- Length: 43.5cm (Total 47.5cm including J-loop handle)
- Diameter: 1.8mm
- EO sterilisation
- Shelf Life: 5 years
- ARTG: 336796

STANDARDS

- ISO 14971:2007 ISO 14971:2007
- ISO11135-1:2007 Sterilisation of health care products-EO

PRODUCT DIAGRAM

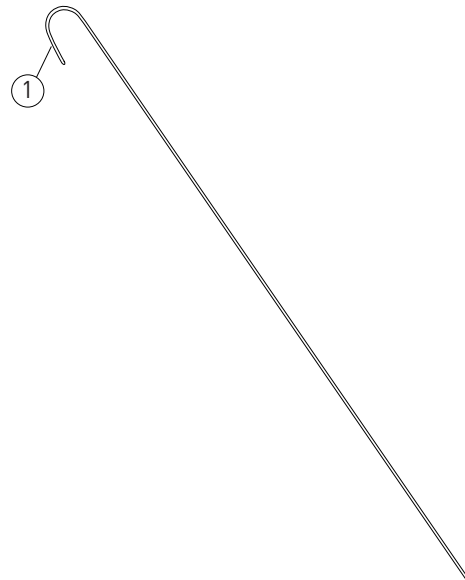

ITEM	DESCRIPTION
1	J-Loop Handle

Cervical Dilating Balloon Introducer

RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
OB020000	Cervical Dilating Balloon Introducer	43.5	20